




NANOPASS® 32.5G 4mm 6mm 8mm

Pen needle – fits all pen-injectors (Type A)



Specifications Nanopass® 32.5G (Ø 0.22 mm)			
Needle length	4 mm	6 mm [▲]	8 mm
Product code	TN*32504M	TN*32506M	TN*32508M
PIP code	381-5321	n.a.	381-5339
Packaging size	100 pieces	100 pieces	100 pieces

▲ Available 2014

"TYPE A"
ISO 11608-2
CE 0197
LATEX
Made of latex-free elastomer

Visit www.terumo-diabetes.com



Terumo Corporation
44-1, 2-chome, Hatagayaya
Shibuya-ku, Tokyo 151-0072,
Japan
Tel.: +81 3 3374 8111
Fax: +81 3 3374 8196

European Headquarters

Terumo Europe N.V.
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 38 12 11
Fax: +32 16 40 02 49

European Offices

**Terumo Europe N.V.
Benelux Sales Division**
Interleuvenlaan 40
3001 Leuven, Belgium
Tel.: +32 16 39 25 80
Fax: +32 16 39 25 99
The Netherlands:
Tel.: 0800 0220396
Fax: 0800 0220414

Laboratoires Terumo France S.A.
Parc Ariane - Bât. Uranus
1, rue Hélène Boucher
78284 Guyancourt Cedex, France
Tel.: +33 130 96 13 00
Fax: +33 130 43 60 85

Terumo Deutschland GmbH
Ludwig-Erhard-Strasse 6
65760 Eschborn, Germany
Tel.: +49 6196 80 230
Fax: +49 6196 80 23 200

Terumo Deutschland GmbH
Zweigniederlassung Spreitenbach
Bodenackerstrasse 3
8957 Spreitenbach, Switzerland
Tel.: +41 56 419 10 10
Fax: +41 56 419 10 11

Terumo Europe España SL
Avda. Juan Carlos I, Nº 13 - 7ª Planta
Edificio Torre La Garena
28806 Alcalá de Henares (Madrid),
Spain
Tel.: +34 9021 01 298
Fax: +34 9021 01 358

Terumo Italia S.r.l.
Via Simone Martini, 143 / 145
00142 Roma RM,
Italy
Tel.: +39 0651 96 14 20
Fax: +39 0650 30 407

Terumo Sweden AB
Nya Varvet, Byggnad 90
Sven Källfelts Gata 18
426 71 Västra Frölunda,
Sweden
Tel.: +46 3174 85 880
Fax: +46 3174 85 890

**Terumo Denmark
Filial of Terumo Sweden AB**
Tel.: +45 7020 93 80
Fax: +45 7020 94 80

Terumo UK Ltd
Tamesis
The Causeway
Egham, Surrey
TW20 9AW, United Kingdom
Tel.: +44 1784 476 200
Fax: +44 1784 476 201

Terumo Russia LLC
BC Northern Tower
13 Floor, 10 Testovskaya Street
Moscow 123317, Russia
Tel.: + 7 495 988 4740
Fax: + 7 495 988 4739

**Terumo Corporation Dubai Branch
(Middle East and Africa)**
Al Masraf Tower, 22nd Floor
P. O. Box 20291 Dubai
U.A.E.
Tel.: + 971 4 221 2220
Fax: + 971 4 221 3330

*Addresses subject to change,
please consult our Website.*





© Registered Trademark
Published by Terumo Europe N.V.
DI 105UK-0813TE-(08.13)E

NANOPASS®

Needle for pen-injectors

Designed to reduce pain

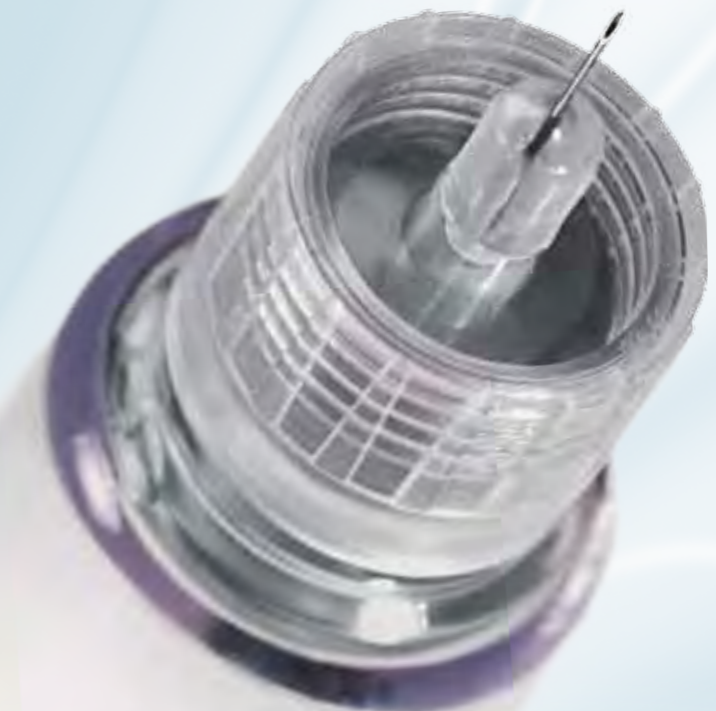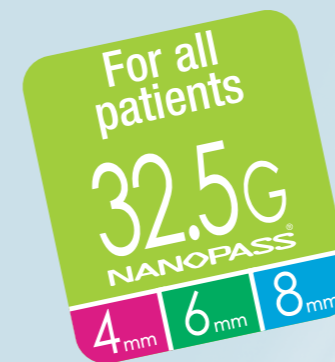

- Unique asymmetric bevel
- Latest technology for improved comfort

Improved quality of life for every patient

NANOPASS®

Less pain, every day

One size made to fit all

- The NANOPASS® 32.5-range offers three different lengths – allowing you to choose the optimal pen needle for each patient type.
- Terumo has used its extensive needle design technology to minimize the pain associated with injections
- Colour coding for easy identification of the preferred pen needle length.



▲ Available 2014

NANOPASS® 32.5 offers patients a more convenient insulin injection by applying dedicated Terumo technologies:


- 1 Asymmetric bevel
- 2 Special thin pen needle size
- 3 Terumo-thin wall technology



The perfect needle for every type of patient

1 Asymmetric bevel


The unique asymmetric bevel of NANOPASS® cuts the skin instead of puncturing it. This decreases the resistance while penetrating the skin[▲] – designed to reduce pain.

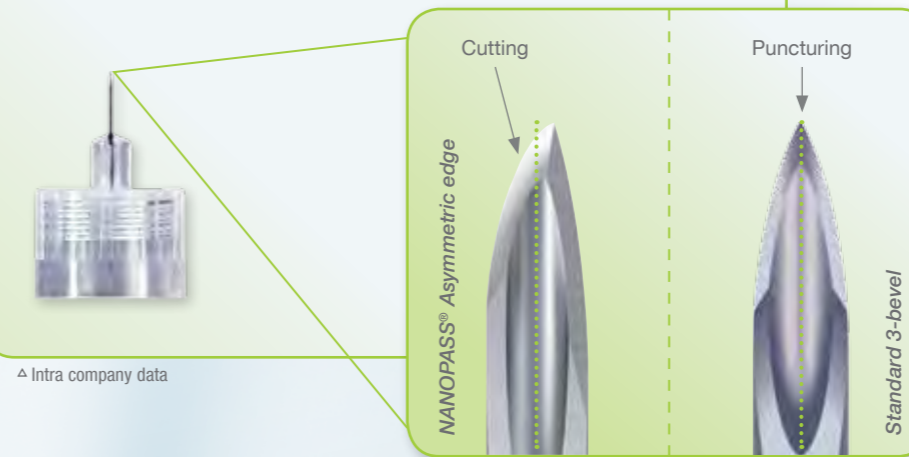

▲ Intra company data

2 Special thin pen needle size

The special thin pen needle size of only 32.5G (0.22 mm) causes less skin trauma than larger needle sizes.

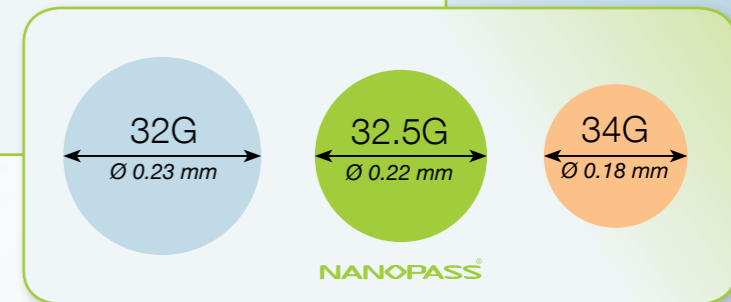

3 Terumo-thin wall technology

The large inner diameter of NANOPASS® 32.5 optimizes the insulin flow for an easy insulin injection.

